Technical Data Sheet

Single-Use Disposable Tourniquets

Description

The tournikidz[®] is a disposable, Single-Use tourniquet designed for assisting venepuncture in blood sampling and intravenous cannulation in patients' under I2 years old.

Function

Strip to be placed around the patients extremity and tapered end to be placed through the pre-cut hole in the opposite (larger) end with the Printed side facing outwards. Adhesive release liner to be removed before tightening the strip to the desired tension to allow for adequate exposure of the desired blood vessels to occur. Strip to be held under tension in situ by sticking the loose end of the strip to the adhesive area exposed underneath. Strip can be easily removed by pulling both ends of the strip away from the adhesive area. Strip should be disposed of after a single use to avoid cross infection risks between subjects.

Product Information

Brand: tournikidz

Manufacturer Part Number: TOURNIKIDZ TSKIDZ4I275DP/ASEPIO0

Product Description: Children's Single Use Tourniquet



Dimensions

50mm Cardboard Core

Materials


Substrate material: TransferPP96 Matt White/SP8000/WG6-M:



Shelf life

Expiry Date (5 years from date of manufacture)



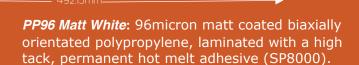

HT236 Adhesive: Tackified water based acrylic adhesive designed to give high tack and adhesion to both polar and non-polar substrates.

Regulatory

Conforms to 93/42/EEC standards as Class I Medical Device

Packaging

200 Individual Strip to be removed from the reel via the self-dispensing box opening by tearing through the pre-perforated cut defining the end of the strip.

Dimensions 95mm x 95mm x 62mm

Materials Printed 500 micron Double sided

white cardboard

