

OG-6 REAGENT, PAP 2A

IVD *In vitro* diagnostic medical device



Cytoplasmic staining reagent acc. to Papanicolaou Counterstain for monochromatic staining of samples in cytology

INSTRUCTIONS FOR USE

REF Catalogue number: 01OG6100 (100 mL) 01OG6500 (500 mL) 01OG61000 (1000 mL) 01OG62500 (2500 mL)

Introduction

OG-6 reagent, Pap 2A is an alcoholic solution of Orange G dye with added phosphotungstic acid (PTA). The first step in using the Papanicolaou staining method implies nuclear staining with a hematoxylin solution, and next two steps consist of contrast staining using the monochromatic OG-6 reagent and one of the polychromatic EA reagent formulations consisting of two acid dyes, the Eosin Y and Light Green SF. The Orange G molecule stains the cytoplasm, and in later stages of the procedure it remains only in the mature, keratinized cells that turn different shades of orange. The third step consists of using one of the polychromatic EA solutions that stains the unstained cellular components, such as squamous cells, nucleoli, cilia, and erythrocytes. Test samples can be gynecological and non-gynecological, such as sputum, urine, and cytological puncture samples. In order to obtain optimal staining results, the properties of OG-6 reagent, Pap 2A are completely in accordance with other Histo-Line reagents used for cytological staining acc. to Papanicolaou - OG-6, Pap 1A reagent, EA 31, Pap 3A reagent, as well as alternative counterstain polychromatic stains, such as EA 50, Pap 3B reagent and EA 65, Pap 3C reagent.

Product description

OG-6 REAGENT, PAP 2A - Counterstain for monochromatic staining of samples in exfoliative cytology. Contains BSC-certified Orange G dye with added phosphotungstic acid and required stabilizers.

Preparing the cytological smear for staining

There are two methods of collecting and preparing the cytological samples:

1. After collecting the cytological sample, place it on the microscope slide, fixate it immediately with a fixative in a spray bottle (CitoSpray), dry it and keep until the staining process. Cytological sample may be fixated and kept until staining by immersing into 95% alcohol solution (Alcohol 95) for a minimum of 30 minutes.
2. Using liquid-based cytology method (LBC) and brush for collecting cytological samples, fixate the sample immediately (CitoFix, CitoFix in transport containers) by removing the brush head and immersing it in the fixative. At the beginning of processing the sample, isolate the cells from the fixative (one of the methods is to centrifuge the fixative) and place them on the microscope slide equally in a single layer. Cytological sample prepared in such a way is ready for staining.

The Papanicolaou staining method, **PROGRESSIVE**

The first stage of staining procedure depends on the method the cytological sample was collected and fixated on the microscope slide.

If the sample is dry and previously fixed using CitoSpray, it is necessary to keep it in a 95% alcohol solution (Alcohol 95) for 10 minutes in order to remove polyglycols. If the section was fixated with a 95% alcohol solution (Alcohol 95), ignore this step. During staining cytology samples (prepared by using the liquid based cytology method (LBC)) that contain low concentration of alcohol, rehydration by descending series of alcohol solutions is not necessary. The procedure starts by rinsing the section using distilled (demi) water and is then stained using Hematoxylin HP, Pap 1A reagent.

1.	Rehydrate in descending series of alcohols (Alcohol 95, Alcohol 80 and Alcohol 70) and in distilled or demineralized water	6-8 dips in each of the 4 exchanges
2.	Staining using Hematoxylin HP, Pap 1A reagent	2-3 minutes
3.	Blue using Scott's solution or Bluing reagent	1 min
	Note: If the mentioned reagents are not available, the section should be blued using indirect stream of water	3-5 minutes
4.	Immerse the sections in distilled/demineralized water.	
5.	Dehydrate in ascending series of alcohols (Alcohol 70, Alcohol 80 and Alcohol 95)	6-8 dips in each of the 3 exchanges
6.	Stain using OG-6, Pap 2A reagent	2-3 minutes
7.	Rinse using 95% alcohol in two exchanges (Alcohol 95)	6-8 dips in each of the 2 exchanges
8.	Stain using EA 31, Pap 3A reagent or EA 50, Pap 3B reagent	2-3 minutes
9.	Rinse using 95% alcohol (Alcohol 95)	6-8 dips
10.	Dehydrate using 100% alcohol (Alcohol 100)	6-8 dips
11.	Dehydrate using 100% alcohol (Alcohol 100)	3-5 minutes
12.	Clear the section in xylene or in a xylene substitute	6-8 dips
13.	Clear the section in xylene or in a xylene substitute	3-5 minutes

Immediately after clearing apply an appropriate HistoMount medium for covering/mounting on the section. If xylene was used, use one of Histo-Line mounting xylene-based media (HistoMount, HistoMount High, HistoMount M, HistoMount DPX, HistoMount C, or universal HistoMount New). If xylene substitute was used, the appropriate covering agent is HistoMount New. Cover the section with cover glass.

Papanicolaou staining method, REGRESSIVE

The regressive staining method creates a better sample differentiation and clearer nuclear structure visibility.

The first stage of staining procedure depends on the method the cytological sample was collected and fixated on the microscope slide.

If the sample is dry and previously fixed using CitoSpray, it is necessary to keep it in a 95% alcohol solution (Alcohol 95) for 10 minutes in order to remove polyglycols. If the section was fixated with a 95% alcohol solution (Alcohol 95), ignore this step. During staining cytology samples (prepared by using the liquid based cytology method (LBC)) that contain low concentration of alcohol, rehydration by descending series of alcohol solutions is not necessary. The procedure starts by rinsing the section using distilled (demi) water and is then stained using Hematoxylin HP, Pap 1A reagent.

1.	Rehydrate in descending series of alcohols (Alcohol 95, Alcohol 80 and Alcohol 70) and in distilled or demineralized water	6-8 dips in each of the 4 exchanges
2.	Staining using Hematoxylin HP, Pap 1A reagent	6 min
3.	Rinse in distilled/demineralized water	6-8 dips
4.	Differentiation using HCL Pap reagent or in 0.1% HCl solution	5-10 seconds
	Note: This step removes excessive hematoxylin from the nucleus and cytoplasm. Discoloration of the nuclei can occur if the section is treated with the differentiation agent for too long.	
5.	Rinse in distilled water	6-8 dips
6.	Blue using Scott's solution or Bluing reagent	1 min
	Note: If the mentioned reagents are not available, the section should be blued using indirect stream of water	
7.	Immerse the sections in distilled/demineralized water.	
8.	Dehydrate in ascending series of alcohols (Alcohol 70, Alcohol 80 and Alcohol 95)	6-8 dips in each of the 3 exchanges
9.	Stain using OG-6, Pap 2A reagent	3 min
10.	Rinse using 95% alcohol in two exchanges (Alcohol 95)	6-8 dips in each of the 2 exchanges
11.	Stain using EA 31, Pap 3A reagent or EA 50, Pap 3B reagent	3 min
12.	Rinse using 95% alcohol (Alcohol 95)	6-8 dips
13.	Dehydrate using 100% alcohol (Alcohol 100)	6-8 dips
14.	Dehydrate using 100% alcohol (Alcohol 100)	3-5 minutes
15.	Clear the section in xylene or in a xylene substitute	6-8 dips
16.	Clear the section in xylene or in a xylene substitute	3-5 minutes

Immediately after clearing apply an appropriate HistoMount medium for covering/mounting on the section. If xylene was used, use one of Histo-Line mounting xylene- based media (HistoMount, HistoMount High, HistoMount M, HistoMount DPX, HistoMount C, or universal HistoMount New). If xylene substitute was used, the appropriate covering agent is HistoMount New. Cover the section with cover glass.

Note

In the case of subsidence in the Hematoxylin HP, Pap 1A solution or formation of metallic glow on the surface, reagent should be filtered before use. Time periods of staining procedures are not completely standardized. The suggested methods are in accordance with Histo-Line reagents' properties and correspond to longtime clinical and laboratory practice. Intensity of staining depends on the period of exposure to stains and reagents. Staining procedure can be changed according to personal preferences if they correspond to the basic principles of cytotechnology.

Results

Blue - nuclei

Yellow-orange - keratinized cells

Pink-red - superficial squamous epithelial cell, erythrocytes, nucleoli, cilia

Green - cytoplasm of all the other cell types (parabasal and intermediate squamous cells, columnar cells, polymorphonuclear leukocytes, lymphocytes, histiocytes, adenocarcinomas, undifferentiated carcinoma cells)

Preparing the sample and diagnostics

Use only appropriate instruments for collecting and preparing the samples. Process the samples with modern technology and mark them clearly. Follow the manufacturer's instructions for handling. In order to avoid mistakes, the staining procedure and diagnostics should only be conducted by authorized and qualified personnel. Use only microscope according to standards of the medical diagnostic laboratory.

Safety at work and environmental protection

Handle the product in accordance with safety at work and environmental protection guidelines. Used solutions and out of date solutions should be disposed of as special waste in accordance with national guidelines. Chemicals used in this procedure could pose danger to human health. Tested tissue specimens are potentially infectious. Necessary safety measures for protecting human health should be taken in accordance with good laboratory practice. Act in accordance with signs and warnings notices printed on the product's label, as well as in Histo-Line material safety data sheet.


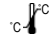



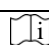


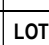

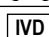


Storing, stability and expiry date

Keep OG-6 reagent, Pap 2A in a tightly closed original package at temperature between 15°C and 25°C. Keep in dry places, do not freeze and avoid exposing to direct sunlight. Date of manufacture and expiry date are printed on the product's label.

References

1. Papanicolaou, G.N. (1941): Some improved methods for staining vaginal smears. J Lab Clin Med.
2. Papanicolaou, G.N. (1942): A new procedure for staining vaginal smears. Science.
3. Carson, F.L., Hladik C. (2009): Histotechnology: A self-instructional text, 3rd ed. ASCP Press.
4. Sherwani, R.K., Khaqn, T. et al. (2007): Conventional Pap Smear and Liquid Based Cytology for Cervical Cancer Screening – A Comparative Study, Journal of Cytology, 24 (4): pp 167-172.

OG6-OT-X, V28-EN12, 31 May 2017, AK/VR

	Refer to the supplied documentation		Storage temperature range		Number of tests in package		Product code		European Conformity
	Refer to supplied instructions		Keep away from heat and sunlight		Valid until		Lot number		Manufacturer
	For in vitro diagnostic use only		Keep in dry place		Caution - fragile				

Histo-Line Laboratories srl
V. G. Di Vittorio 30
20090 Pantigliate (MI)
ITALY
www.histoline.com

