



# HIV Rapid Test

Screening Test for HIV Infection



# HIV 1/2/O Rapid Test

HIV (human immunodeficiency virus) is a lentivirus that causes HIV infection and over time acquired immunodeficiency syndrome (AIDS). Infection with HIV occurs by the transfer of blood, semen, vaginal fluid, or pre-ejaculate. It may also be transmitted from mother to fetus during pregnancy or through breast feeding.

#### **General Information**

ACON HIV 1/2/O Rapid Test is a rapid chromatographic immunoassay for the qualitative detection of antibodies to HIV-1, including subtype O, and HIV-2 simultaneously in serum, plasma or whole blood.

- Detection: IgG & IgM (by recombinant HIV-1 (gp41), HIV-2 (gp36) and HIV-O Ag)
- Specimen: Serum/Plasma, Serum/Plasma/Whole Blood
- Sensitivity: >99.9%
- Specificity: >99.0%
- Time to Result: 10-20 minutes
- Shelf life and storage temperature: 24 months at 2-30°C

#### **Materials Provided**

- Test Strip / Cassette
- Dropper
- Buffer (for S/P/WB test only)
- · Package insert

## **Test Procedure and Interpretation**

## **HIV Ab Strip**


- For S/P test: Add 2 drops (80μL) of specimen
- For S/P/WB test: Add 2 drops (50 $\mu$ L) of specimen, then add 1 drop (40 $\mu$ L) of buffer
- Read result at between 10-20 minutes

#### **HIV Ab Cassette**


- For S/P test: Add 2 drops (80µL) of specimen
- For S/P/WB test: Add 2 drops  $(50\mu L)$  of specimen, then add 1 drop  $(40\mu L)$  of buffer
- Read result at between 10-20 minutes

#### HIV Ab Tri-line Cassette

- Add 2 drops (50µL) of specimen, then add 1 drop (40µL) of buffer
- Read results at between 10-20 minutes

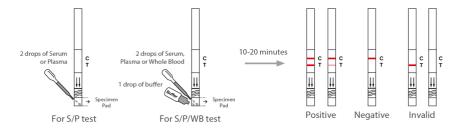



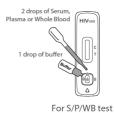







2 drops of Serum, Plasma or Whole Blood HIV







10-20 minutes



## Ordering Information

Product Name	Catalogue	Package
HIV 1/2/O Rapid Test Strip (Serum/Plasma)	L031-10111	50 Tests/Kit
HIV 1/2/O Rapid Test Strip (Serum/Plasma/Whole Blood)	L031-10121	50 Tests/Kit
HIV 1/2/O Rapid Test Cassette (Serum/Plasma)	L031-10131	40 Tests/Kit
HIV 1/2/O Rapid Test Cassette (Serum/Plasma/Whole Blood)	L031-10141	40 Tests/Kit
HIV 1/2/O Tri-line Rapid Test Cassette (Serum/Plasma/Whole blood)	L031-11411	40 Tests/Kit



www.aconlabs.com

ACON Laboratories, Inc.

5850 Oberlin Drive, #340, San Diego, CA 92121, USA

Tel: 1.858.875.8000 Fax: 1.858.200.0729 Email: info@aconlabs.com

1150971104 © 2023 ACON Laboratories, Inc.