



Quantitative Real-time PCR system

.....

(ZIP-96V)



Molecular




Quantitative Real-time PCR system (ZIP-96V)

Quantitative real-time PCR system ZIP-96V is 96-well real-time PCR instrument based on Fresnel lens optical signal acquisition technology, time-resolved signal separation technology and unique temperature control technology. It is designed for multiple fields of basic research in medicine and biology, such as pathogen detection, genotyping, gene therapy drug, gene expression, public health, etc.




Features

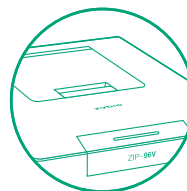


User-friendly for first-time and experienced operators

- No ROX correction required
- Preset projects, such as COVID-19, influenza A/B project

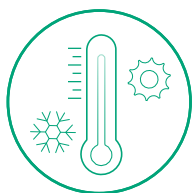

Advanced temperature control block

- The high-level peltier chip combined with a new generation of semiconductor heating and cooling technology greatly improves the heating and cooling rate, and the heating rate can reach up to 6°C/s
- Six independent temperature control blocks improve the accuracy and uniformity of temperature control

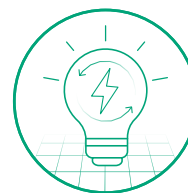

Excellent optical system


- Use high efficiency LEDs light source, maintenance-free for life
- Unique test tube side daylighting technology, four-channel fast daylighting only takes 20s, automatically adjusts gain, improves the fluorescence signal sensitivity and signal to noise ratio
- Unique time-resolved signal separation technology, no multi-color crosstalk

Multiple Applications

- Basic scientific research
 - Gene expression
- Pathogen detection
 - Public health
- Genotyping
- Gene therapy drug

Four fluorescence Channels, support most of common dyes

Detection Reporter dye:

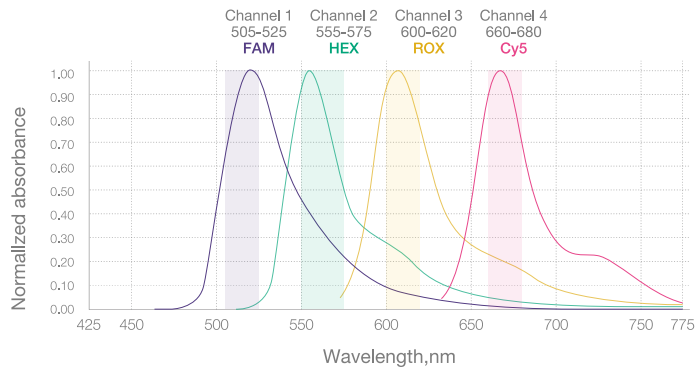

Technical parameters

| Basic parameters | | Optical system | |
|---------------------------|---------------------------------|----------------------------|---|
| Model | ZIP-96V | Excitation light source | High efficiency LEDs |
| Dimensions(W×D×H, mm) | 425*320*205 | Detector | Photodiode |
| Weight | 15kg | Excitation wavelength | 460nm~480nm, 525nm~545nm, 575nm~595nm , 600nm~620nm |
| Throughput | Up to 96 tests per run | Emission wavelength | 505nm~525nm, 555nm~575nm, 600nm~620nm , 660nm~680nm |
| Reaction volume | 20-120μL | Fluorescence repeatability | < 2% |
| Sample Format | 0.2ml PCR tube/ | Fluorescence linear | r > 0.99 |
| | 8-strip PCR tube/ 96-well plate | | |
| | Channel 1: FAM | | |
| | Channel 2: VIC, HEX | | |
| Supported dye | Channel 3: ROX | | |
| | Channel 4: CY5 | | |
| Fluorescent channels | 4 detection channels | | |
| Thermal cycle technology | | Software | |
| Temperature control | Peltier | Applications | Qualitative/ Absolute Quantification Analysis |
| Temperature control block | 6 blocks | Scan mode | All-channel scanning |
| Temperature range | 4~99℃ | Data export | Excel |
| Maximum heating rate | ≥2.5℃/s | Report | Multiple printing templates/ customized template |
| Maximum cooling rate | ≥2.0℃/s | | |
| Temperature fluctuation | ± 1℃ | | |
| Temperature accuracy | ≤ ±0.5℃ | | |
| Hot-lid temperature range | 30~110℃ | | |
| | | Working Environment | |
| | | Power supply | 100-240 V AC, 50-60 Hz |
| | | Operating system | Windows 7/10 |
| | | Interface | USB |
| | | Working Temperature | 10-30℃ |
| | | Working Humidity | 20-85 % |

Classic Examples

SARS-CoV-2 Nucleic Acid Detection (PCR-Fluorescent Probe Method) qPCR technology was used to detect and diagnose suspected case of COVID-19.

- **Detection kit:**

SARS-CoV-2 nucleic acid detection kit(Zybio)

- **Isolation kit and system:**

nucleic acid extraction kit(Zybio), nucleic acid isolation system EXM3000(Zybio)

- **Assay method :**

Take 200μL of sample(throat swab preservation solution) to be tested, and extract RNA from the sample with nucleic acid extraction kit and isolation system EXM3000. According to the instruction of nucleic acid detection kit, mix extracted RNA, the positive control product, and the negative control product with PCR reaction solution and enzyme respectively, mix well and centrifuge and then do amplification and detection.

- **Amplification parameters :**

| Steps | | Temperature | Time | Cycle |
|-------|--|-------------|--------|-------|
| 1 | UNG reaction | 37°C | 1 min | 1 |
| 2 | Reverse transcription | 50°C | 5 min | 1 |
| 3 | Initial denaturation | 95°C | 2 min | 1 |
| 4 | Denaturation | 95°C | 5 sec | 45 |
| 5 | Amplification and fluorescence detection | 60°C | 30 sec | |

Fluorescence Detection: Step 5.

Report Fluorescence Setting: FAM, ROX, VIC;


Quenching Fluorescence Setting: None;Passive Reference Setting: None.

- **Cut-off Value:**

result is claimed positive when the Ct of two targets(FAM, ROX) <40.

- **Result analysis:**

Positive control and negative control are normal, sample to be tested is positive.The corresponding amplification curves are as follows :

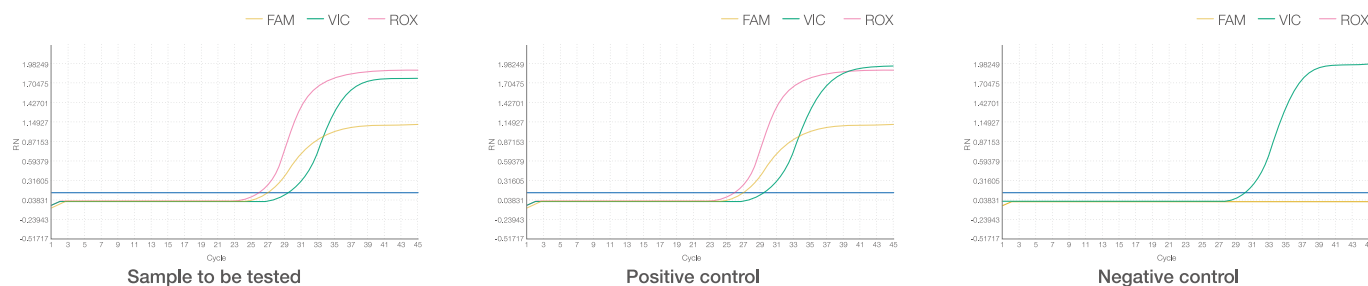



Zybio Inc.

Address: Floor 1 to Floor 5, Building 30, No.6 of Taikang Road,
Block C of Jiannqiao Industrial Park, Dadukou District, Chongqing, China 400082
Tel: +86-23 6865 5509 Fax: +86-23 6869 9779
Email: info@zybio.com Website: www.zybio.com

EN-C-FZ-ZIP-96V-20220414H

Molecular