

Automated Enzyme Immunoassay System

# AIA-360

## SPECIFICATIONS

<b>Assay Principle</b>	Fluorescence Enzyme Immunoassay
<b>Detection</b>	LED illuminant, non flow cell / TOP-TOP photometry method
<b>Throughput</b>	36 tests per hour
<b>Reaction time</b>	Antigen antibody reaction: 10 minutes
<b>Sample loading Capacity</b>	Maximum 25 samples
<b>Reagent loading Capacity</b>	Maximum 25 tests
<b>Test for each sample</b>	4 tests per sample
<b>Sampling</b>	Fixed probe with clot detection function
<b>Sample containers</b>	Primary tube : 13 × 75, 13 × 100, 16 × 75, 16 × 100 mm (Diameter × Length)  Dedicated sample cup
<b>Specimen barcodes</b>	CODE39, CODE128, ITF, NW-7 compatible (compliant with ASTM14466-92 and CLA guideline compliant barcodes) Min. element width (narrow) : 0.191 (0.254 or over recommended)
<b>Interface</b>	RS-232C
<b>Size / Weight</b>	400 (W) × 400 (D) × 520 (H) mm / 29 kg
<b>Power requirement</b>	AC100-240 V 50 / 60 Hz 250 VA
<b>Operating Environment</b>	Temperature: 15-30 °C, Humidity: 40-80 %
<b>Conformance Status</b>	
EMC standard	EN60601-1-2: 2001
CISPR class and group categories	Class A, Group 1
FCC	Part15 Sub, Part B Class A

AIA-PACK, AIA are registered trademarks of Tosoh Corporation.



TOSOH

\* Specifications and design may be changed without prior notice.

## TOSOH CORPORATION Bioscience Division

3-8-2 Shiba, Minato-ku, Tokyo 105-8623, Japan  
Phone : +81-3-5427-5181  
Fax : +81-3-5427-5220



# Biosite

MEDICAL INSTRUMENTS

MANILA Cebu Davao  
Sales Representative: 09171115008  
Technical Services Department: 09171115008

Printed in Japan 5209GX[Cat.No.0992219]D



TOSOH

Automated Enzyme Immunoassay System

# AIA-360

Compact • Fast • Simple • Accurate  
Immunoassay near-at-hand



# AIA-360

# Immunoassay near-at-hand



## Compact and Lightweight

The AIA-360 requires little counter space. Its compact design is suitable for physician office laboratories, work stations for small hospitals, cardiac and specialty testing.



Depth 400 mm  
Height 520 mm  
Width 400 mm

## Easy to Use, Easy to Load

Once samples and test cups are loaded, the operator presses START and walks away.




The operator may load bar-coded primary tubes or sample cups directly on the system. Up to four assays per specimen can be loaded.

The touch screen user interface provides easy access to screens associated with daily operations.



## Fast reporting

With a throughput of 36 tests per hour, the AIA-360 reports its first result in approximately 19 min.

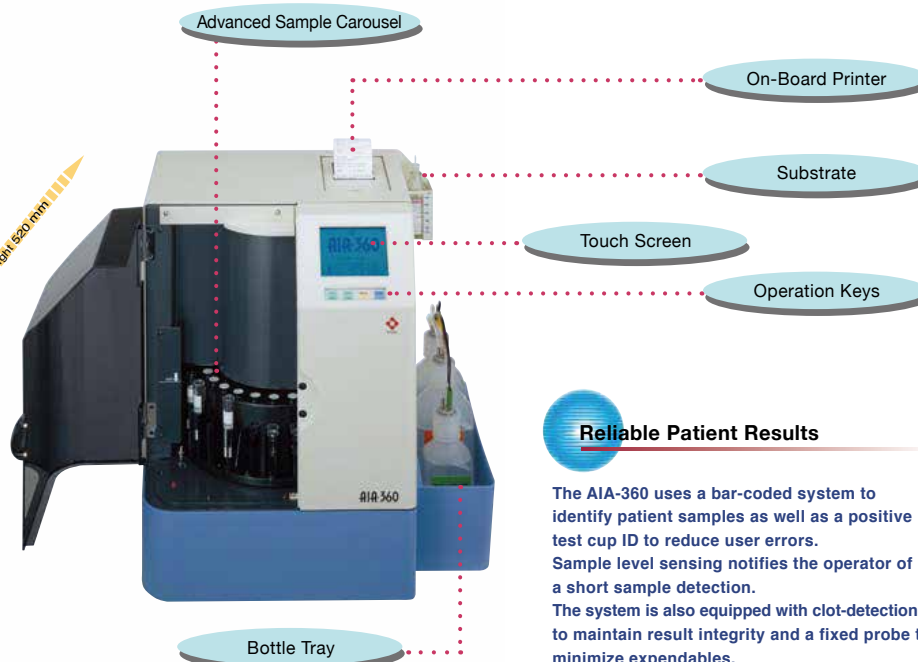

## Continuous Random Access

The Sample carousel holds up to 25 patient samples and 25 test cups at once. The AIA-360 allows the user the option loading the entire carousel or adding the samples as needed. Additional samples and tests can be added at any time without interrupting sample processing.



## Reliable Patient Results

The AIA-360 uses a bar-coded system to identify patient samples as well as a positive test cup ID to reduce user errors. Sample level sensing notifies the operator of a short sample detection. The system is also equipped with clot-detection to maintain result integrity and a fixed probe to minimize expendables.

## Wide array of Parameters

With Tosoh's expanding test menu including cardiac markers, diabetes markers, tumor markers, thyroid and reproductive hormones and more, the AIA-360 can run many test menu for ST AIA-PACK.



Automated Enzyme Immunoassay System


# AIA-360

## Simple Assay Operation

### OPENING Screen



### ASSAY MONITOR Screen

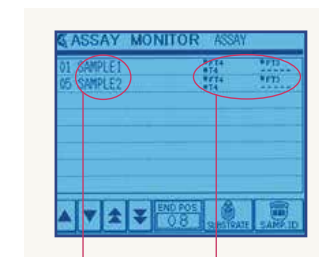

**SAMPLE ID**  
Sample ID (up to 16 characters) can be displayed on the screen when the specimen is read.

**ANALYTE**  
Analytes in the reagent cups read from the camera are displayed for each sample.

### Sample Assay Result Printout

