

Description

Moulded dispenser to hold 150 tournistrip tourniquets.


Function

The dispenser can be wall mounted, stuck to a flat surface using an adhesive or fixed to a medical grade trolley. When refilling is required the dispenser can be accessed by pressing the latch on the front lid moulding (see diagram). This will unlock the lid and enable the front lid to open. The refill pack can then be inserted and the 1st strip of the pack to be pushed through the aperture.

Dimensions—138 mm x 140 mm

Materials


Polypropylene

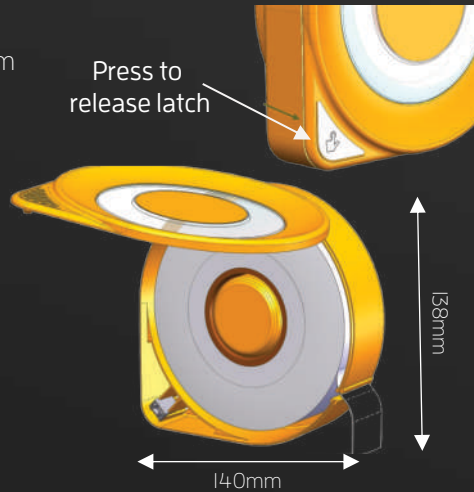

Brand: tournistrip Dispenser with Refill

Manufacturer Part Number:
508OLS692D/ASEPI50

Product Description: Tournistrip: Refill Pack of 150 Tournistrisps

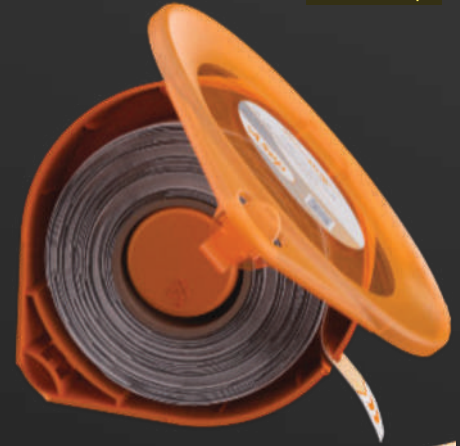


Description

150 Tournistrip eco® refill for the tournistrip dispenser.

Function

Users unwrap the liner and place the refill inside the dispenser on the centrally located holder. Unwind the 1st tourniquet on the strip and place through the aperture of the dispenser.

Dimensions

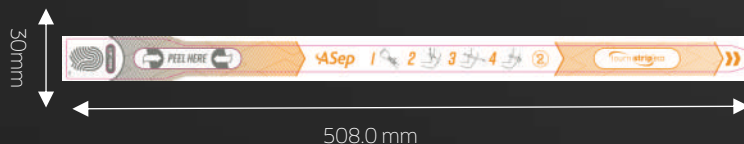

Materials

Base Material: A white bi-axially oriented, polypropylene film with a printreceptive topcoat. Basis Weight 5 6 g/m2 I S O 536 Caliper 57 pm ISO 534

Adhesive: A general purpose S692N permanent, acrylic based adhesive.

Liner: BG40 white, a supercalendered glassine paper. Basis Weight 61 g/m2 I S O 536 Caliper 5 4 pm I S O 534

Laminate: Total Caliper 130 pm±10% I S O 534

Shelf Life

Expiry Date (5 years from date of manufacture) Displayed on packaging

Product Information

Brand: tournistrip Dispenser Refill Pack

Manufacturer Product Code: 508OLS692/ASEPI50

Product Description: Tournistrip Dispenser Mount with 150 Tournistrip eco™ Refill

