

Making time for smiles



Moveable phlebotomy assistant




Specifications

Key information	Specifications	Key information	Specifications
Model name	HENc	Available tube	13/75, 13/100, 16/100
Dimensions	515(W) x 970(D) x 1077(H) mm	Max. variation of tube type	8 different types of tube (50 tubes per type)
Components	HENc, Control software	Tube capacity	400 ea in total
Power source	AC 100-240V 50/60Hz	Compatible tube brand	BD, GREINER, AMPULAB, INSEPACK, IMPROVE
Interface	TCP/IP, RS232, Wi-Fi(optional)	Label capacity	Thermal label 2000 labels per roll
GUI Display	3.5" Touch display	Printing method	Thermal printing
Available barcodes	Code39, Code93, Code128, Codabar, ITF	International standard	KC
Tube supply method	Direct tube pallet supply		

*Above specifications are subject to change without notice due to manufacturer's reason.

Dimensions

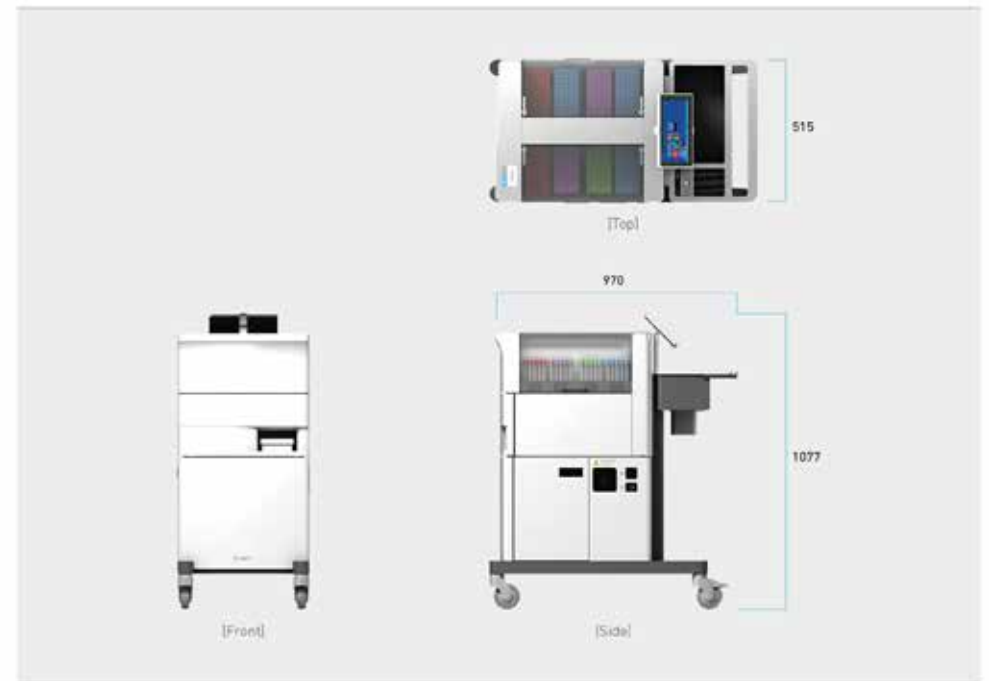

Contact us

ENERGIUM

3ba-610 Sihwa Industrial Complex, 2189-10 Jeongwang-dong, Siheung-si, Gyeonggi-do, South Korea

Tel. +82 2 6919 0311 Fax. +82 2 6919 0300 E-mail, sales@energium.com

www.energium.com



Biosite
MEDICAL INSTRUMENTS

MANILA 201 C&B Bldg. 2nd Floor, San Juan City 1505, Metro Manila, Philippines
CEBU 101 E. 1st St. 2nd Floor, Maguikang, Cebu City 6000, Philippines
DAVAO 101 E. 1st St. 2nd Floor, Maguikang, Davao City 8000, Philippines
 E-mail Address: customerservice@biositeph.com
 Customer Care Hotline: 09171115008
 Sales Department: 09171115008
 Technical Service Department: 09171114008
 @customerservicebiositeph @biositemedicalinstruments @biositeph.com @biositemedicalinstruments

Moveable



HENC is mobile, allowing it to be taken out of the phlebotomy lab and used in regular wards as well.

About HENC

HENC automatically provides medical technicians with the essential items required to perform a phlebotomy (tube with specimen label attached, specimen labels, patient identification labels) based on the phlebotomy information recorded in the Hospital Information System (HIS, LIS, OCS, EMR)



Tablet PC



A PC with the HENC program installed helps medical staff carry out phlebotomy and record the results without error.

HENC Program



A program for easy and convenient checking and recording (information update) of phlebotomy information.



Roller mixer



Assists in uniformly mixing samples in tubes with additives.

Needle disposal container tray



Needle disposal container can be attached to and used with HENC.

Cabinet



Phlebotomy tools and spare tube pallets can be stored.

HENC package

Added options

HENC benefits



HENC



HENC Software



Barcode scanner



Prevent medical errors



Cut waiting times



Cost reduction



Improved convenience



Improved image



Biosite
MEDICAL INSTRUMENTS

MANILA
305 Old Bonyon Avenue,
San Juan City 1503
Metro Manila, Philippines
customer.service@biositeph.com

Customer Care Hotline:
Technical Services Department: 09171114008
Sales Department: 09171115008
Email Address:
customer.service@biositeph.com